

TRAF3IP2 Rabbit Monoclonal Antibody

OPM18454YP

Source	Rabbit	Storage	-20°C, Avoid freeze / thaw cycles
Isotype	IgG	Applications	WB, ELISA


Note: Centrifuge before opening to ensure complete recovery of vial contents.

Product Information

Buffer	PBS with 0.02% sodium azide and 50% glycerol, pH 7.3
Dilution	WB 1:1000-1:4000
Purification method	Protein A purification
Observed MW	60-70kDa

Other applications have not been tested. Optimal dilutions/concentrations should be determined by the end user.

Image

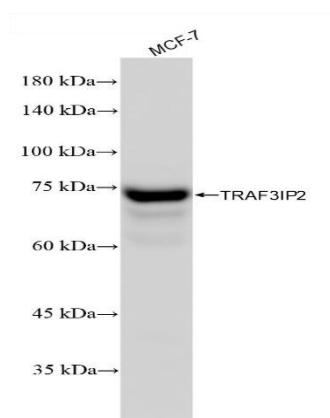

MCF-7 cells were subjected to SDS PAGE followed by western blot with TRAF3IP2 Antibody at dilution of 1:2000 incubated at room temperature for 1.5 hours.

Research Use

For research use only, not for use in diagnostic procedures

Legend

Applications: WB-Western Blot; IHC-Immunohistochemistry; IF-Immunofluorescence; IP-Immunoprecipitation; FC-Flow cytometry; ChIP-Chromatin Immunoprecipitation

Reactivity: H-Human; R-Rat; M-Mouse; Mk-Monkey; Dg-Dog; Ch-Chicken; Hm-Hamster; Rb-Rabbit; Sh- Sheep; Pg-Pig; Z-Zebrafish; X-Xenopus; C-Cow.

Please contact Origin Diagnostics and Research for further assistance

www.originlab.in

info@originlab.in

+91-7736237778