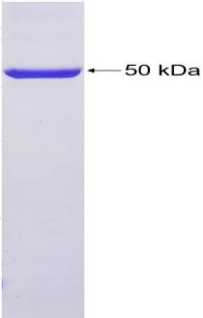


Recombinant Human COX2 protein His-tagged

COX2-2698H

Product Overview	Recombinant Human COX2 protein(20-366 aa), fused with N-terminal His tag, was expressed in E. coli.
Species	Human
Source	E.coli
Protein Length	20-366 aa
Tag	N-His
Form	The purified protein was Lyophilized from sterile PBS (58mM Na ₂ HPO ₄ ,17mM NaH ₂ PO ₄ , 68mM NaCl, pH7.4). 5 % trehalose and 5 % mannitol are added as protectant before lyophilization.
Purity	85%, by SDS-PAGE with Coomassie Brilliant Blue staining.
Storage	Short-term storage: Store at 2-8°C for (1-2 weeks). Long-term storage: Aliquot and store at -20°C to -80°C for up to 3 months, reconstitution with sterile water and addition of an equal volume of glycerol. Avoid repeat freeze-thaw cycles.
Reconstitution	Reconstitute at 0.25µg/µl in 200µl sterile water for short-term storage. After reconstitution with sterile water, if glycerol has no effect on subsequent experiments, it is recommended to add an equal volume of glycerol for long-term storage.
AA Sequence	PCCSHPCQNRGVCMSVGFQYKCDCTRTGFYGENCSTPEFLTRIKLFLKPTPNTVH YILTHFKGFWNVVNNIPFLRNAIMSYVLTSRSHLIDSPPTYNADYGYSWEAFSNLSY YTRALPPVPDDCPTPLGVKGKKQLPDSNEIVEKLLLRKFIQDPQGSNMMFAFFAQH FTHQFFKTDHKRGPAFTNGLGHGVLDLNIHYGETLARQRKLRLFKDGMKYQIIDGE M YPPTVKDTQAEMIPPQVPEHLRFAVGQEVFGLVPGLMMYATIWLREHNRVCDVLK QEHPEWGDEQLFQTSRLILIGETIKIVIEDYVQHLSGYHFKLKFDPELLFNKQFQYQN RIAAE
SDS-PAGE	 <p>← 50 kDa</p>

GENE INFORMATION

Gene Name	MT-CO2 mitochondrially encoded cytochrome c oxidase II [Homo sapiens (human)]
Official Symbol	MT-CO2
Synonyms	COX2; MT-CO2; COII; COXII; MTCO2
Gene ID	4513
UniProt ID	P00403